

Straumann® n!ce®

Fully crystallized glass-ceramic.
n!ce® to meet you.




lithium
Li
6.941
silicium
Si
28.085
oxygen
O
15.999



Constantly innovating to improve patient outcomes.

Through continuous innovation, Straumann® has redefined the boundaries of tooth replacement for dental professionals and their patients.

Our newest innovation, Straumann® n!ce®, is a proprietary glass-ceramic material that combines biocompatible and esthetic materials to create novel, fully crystallized milling blocks for dental professionals.

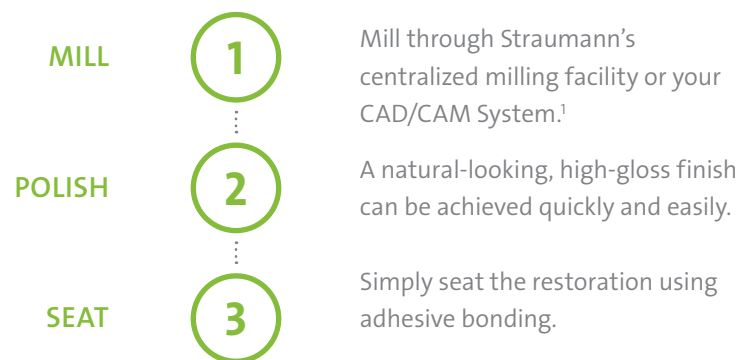

The result is a highly esthetic, natural-looking restoration.

¹ CEREC® Milling or Straumann M Series (coming soon)



Meet and exceed patient expectations

Not only are you helping to build natural-looking, confident smiles with Straumann n!ce, but you can feel confident in Straumann's reliability and quality.

Straumann n!ce allows you to customize each restoration to meet the needs of your patient. Whether you are looking to mill in house or send to Straumann's centralized milling facility, Straumann n!ce is the perfect solution for milling crowns (including hybrid abutment crowns), partial crowns, inlays, onlays and veneers.

Customize every smile:

- Available in HT and LT
- Bleach, A1, A2, A3, B2, B4, and C2 shades are available, in correspondence with the VITA classic shade guide

n!ce restorations can be seated using adhesive or self-adhesive cementation. Polish, stain and glaze as necessary.



¹ Use for exempt applications, specifically single tooth and single abutment (e.g. one-piece stock or CARES® Abutments) borne restorations only. This does not include direct-to-implant applications or use with Ti bases.

Benefit from cost and time savings

n!ce® is a lithium aluminosilicate ceramic reinforced with lithium disilicate. It is a fully crystallized material that does not require any post-mill sintering. It can be ground, polished, and seated immediately, saving time and money.

PROPERTIES

Description	Value
Flexural strength [MPa]*	350 ± 50
CTE (100 – 500 °C) [10 ⁻⁶ K ⁻¹]*	7.1 ± 0.5

* According to ISO 6872

UNIQUE PROPRIETARY
GLASS-CERAMIC DEVELOPED
AND MANUFACTURED
BY STRAUMANN®

A natural surface finish
It's n!ce[®] to mill a smile

NATURAL APPEARANCE
WITH STRENGTH
AND LONGEVITY.

Straumann n!ce requires no additional investment in hardware or software, as the material can be easily processed using standard milling equipment.¹

The unique consistency of the materials allows for simpler milling and handling. With its smooth surface finish and margin lines, the milled element can be polished to achieve a high-gloss restoration. The high-quality, fully crystallized material guarantees a precise fit combined with strength that is easy to use.

From milling to n!ce smiles
in just a few steps.

; -)



¹ C14 blocks for CEREC® Milling or Straumann M Series (coming soon)

Straumann North American Headquarters

Straumann USA, LLC

60 Minuteman Road

Andover, MA 01810

Phone 800/448 8168 (US) • 800/363 4024 (CA)

Fax 978/747 2490

www.straumann.us • www.straumann.ca