

# SAVE TIME - SAVE MONEY - REVOLUTIONIZE YOUR PRACTICE

Straumann brings together the entire dental team with our full, validated digital process designed to simplify the practice of implant dentistry from start to finish.

## STRAUMANN IMPLANTS. STRAUMANN ABUTMENTS. STRAUMANN AUTHENTICITY.

WOULD YOU SETTLE FOR ANYTHING LESS?




1. Intra-Oral Scan of Patient  
Example: iTero®



2. Order iTero Models  
3. Send Model to Lab



4. gonyX™ Laboratory Fabricates Scan Appliance

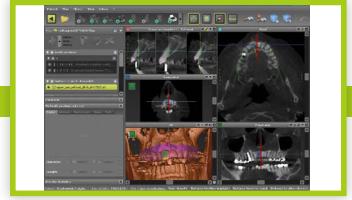


5a. Patient Scanned with Scan Appliance – DICOM Data Transferred to coDiagnostiX™ Surgical Planning Software

5b. iTero Scan Imported to coDiagnostiX™ - Allows Clinician to See Bone and Tissue



6a. Surgical Planning  
6b. Lab Fabricates Surgical Guide Based on Surgical Plan

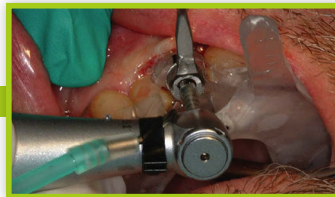

7a. Guided Surgery – Placement of Implant

7b. Second iTero Scan with Mono Scanbody in Straumann® Implant



8. Order iTero Milled Model with Repositionable Implant Analog



9. Import Mono Scanbody STL Data into CARES® Visual Software to Virtually Design Customized Abutment and Temp/Final Restoration

Images courtesy of William L. Wiley, DDS and Champions Crown and Bridge